

# KerrHawe OptiClean™

Finally an easy and efficient way to remove temporary cement.

- ▶ **Faster, easier and complete removal of temporary cement.**
- ▶ **Excellent access due to its specific design.**
- ▶ **Clear view of the working surface since no paste or pumice slurry is required.**
- ▶ **Gentle to the soft tissue.**



## ▶ Colour coding:

KerrHawe offers comprehensive colour systems. OptiClean follows the colour coding of OptiDisc, the first translucent finishing and polishing disc (Light blue: "fine" abrasivity).



## Order information

Art. No.	Type of Package	Contents:
950	OptiClean Kit	12 OptiClean instruments


P11900 E-09

KerrHawe SA  
P.O. Box 268  
6934 Bioggio  
Switzerland

Freephone: 00800-41-050-505  
Fax: +41-91-610-0514  
[www.KerrHawe.com](http://www.KerrHawe.com)



Dentists have been waiting for this for ages.

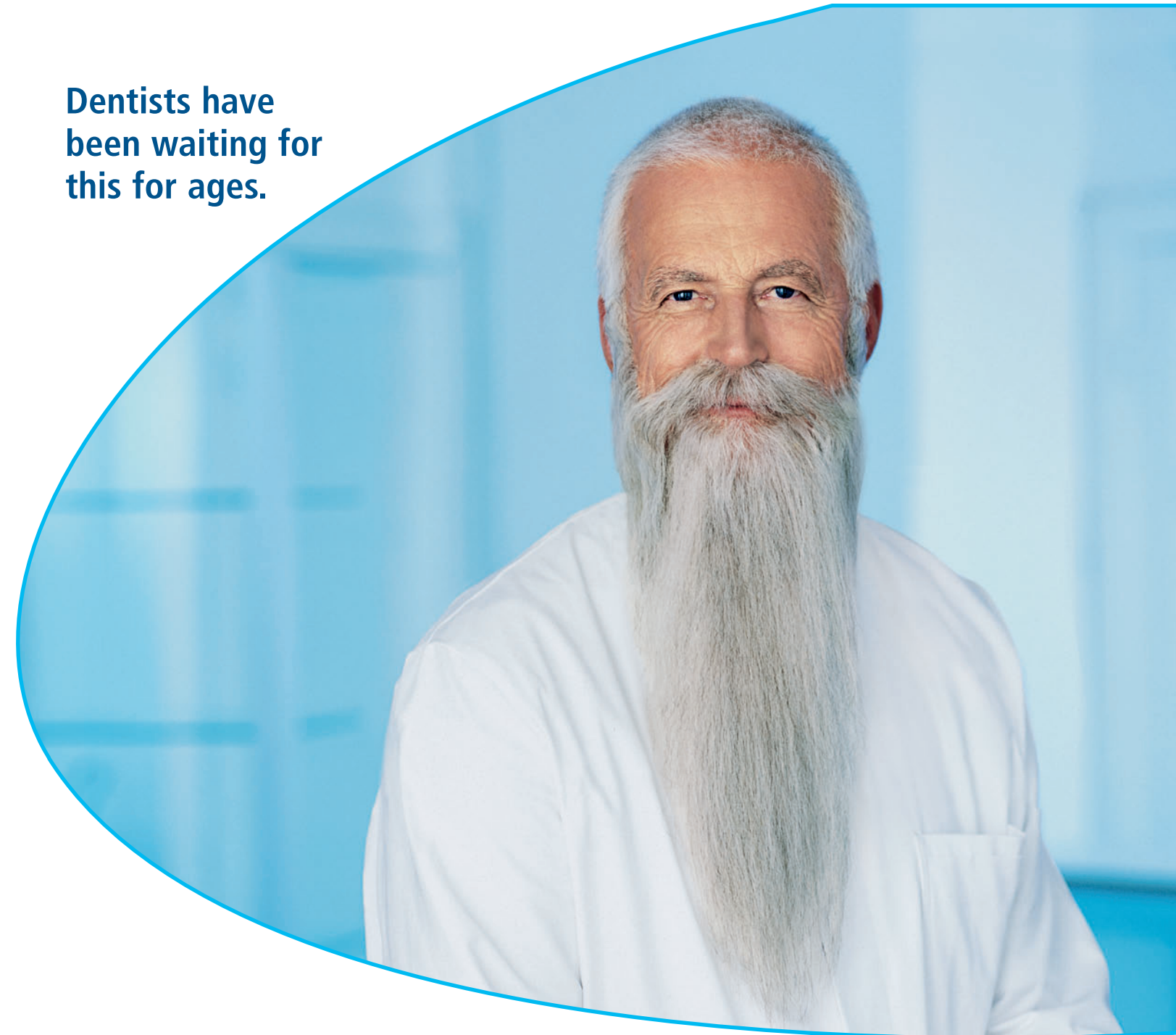

# KerrHawe OptiClean™

The first rotary instrument for the removal of temporary cement and debris.



New in the KerrHawe system for temporary restoration:



## Our innovation for optimal dental care: The first rotary instrument for the removal of temporary cement and debris.

### Great features:

#### Easy and efficient handling:

- ▶ Faster, easier and complete removal of temporary cement.
- ▶ Excellent access due to its specific design.
- ▶ Clear view of the working surface since no paste or pumice slurry is required.

#### More gentle than common procedures:

- ▶ Low risk of abrasion on the preparation and adjacent teeth.
- ▶ Gentle to the soft tissue.
- ▶ Highly hygienic because it's single use only.



US Pat. 5 882 201

So far there has been no efficient way to remove temporary cement. But now you have KerrHawe OptiClean: indicated for the removal of temporary cement and debris.



▶ Removal of temporary restoration from preparation.



▶ Temporary cement removal with KerrHawe OptiClean.




▶ Perfect result after using KerrHawe OptiClean.

### KerrHawe OptiClean technical details:


- ▶ Optimised abrasive particles for efficient cement removal and minimal abrasion on dentine.
- ▶ A-silicone matrix for excellent flexibility and dynamic stress behaviour.
- ▶ High-performance aromatic polymer guarantees perfect and secure fit of the mandrel in the handpiece.
- ▶ Smallest tool available on the market with flexible coating and 1.65 mm tip diameter.



2 mm  
Magnification: 17x

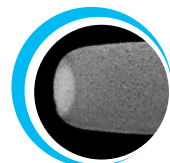

500 µm  
Magnification: 55x

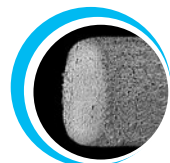

200 µm  
Magnification: 150x

Field Emission SEM Study by Prof. P. Lambrechts, Catholic University of Leuven, Belgium



**George Salem, DMD, FAGD**  
Braintree, Massachusetts, USA

“OptiClean is perfect in its execution! It removes all traces of temporary cement and delivers a perfectly clean cementation surface of the *entire* tooth preparation in literally seconds. I love using OptiClean because it is quick, easy, provides a higher quality surface for cementation and therefore provides a higher quality restoration for my patients.”

Clinical photos by Dr. med. dent. M. Brandenbusch, Johannes Gutenberg University Mainz, Germany