

BRILLIANCE

what you see matters



Eighteenth Brilliance series Surgical Loupes combine Schott Glass and Achromatic Optical Path Design to assure high image quality. The Flip-Up function can bring unobstructed vision and eye contact with your patients with a single adjustment. The whole structure made of ultra-light alloy creates a cool metallic look, ensuring excellent structural strength and durability, while minimizing the weight and the pressure on the nose when wearing.

Intelligent damping torque effect
Easy to adjust and with excellent stability

Super light frame
Ergonomic design with optional materials and colors


Protective shield
Besides protective shields, prescription lenses can also be fitted by an optometrist

Schott Glass
Ensure high resolution, clear and sharp field of view

Refined lens coatings
Multi-layer coating technology, Anti-fog, Anti-scratch, Anti-glare, to ensure more visible light through the glass

Asynchronous PD adjustment
Laser printed scale lines, accurate and intuitive

Aerospace-grade aluminium magnesium alloy
Our loupes are exceptionally high strength and light

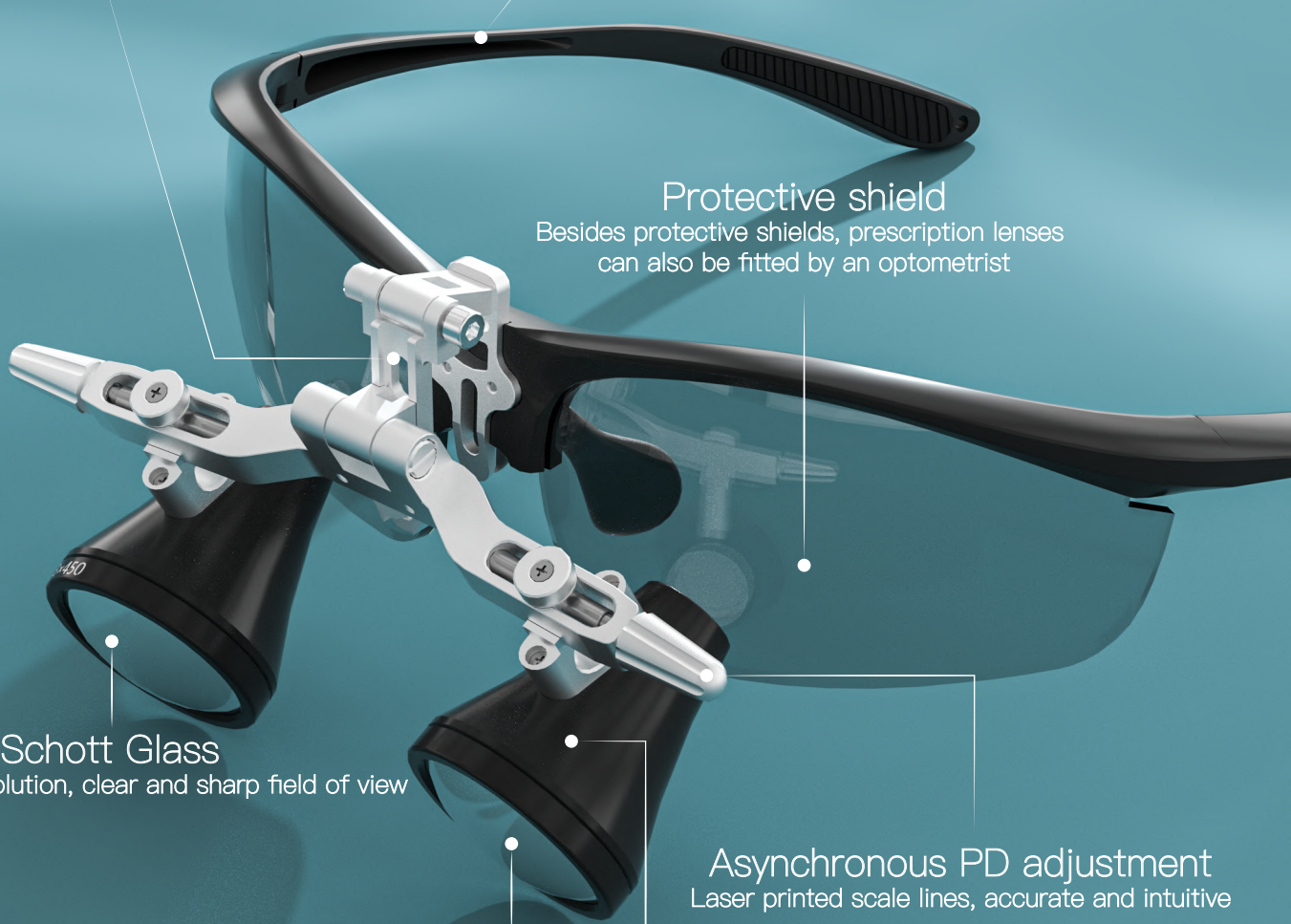

Two Options!



Two hinges type



Three hinges type




Both can be quickly and accurately reset after flipped-up without extra adjustment. Get into the next use immediately, and greatly improve the work efficiency.

Your ideal first pair of Loupes!

The 3x magnification loupes offer extreme versatility and are perfectly suited for general dentistry and those early on in their careers.



Frame Material: AlMg, TR90

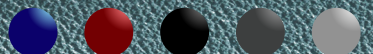Frame Color: 

Magnification	Working Distance	Field of View	Depth of Field	Net Weight
3.0x	350mm	68mm	≥90mm	71g
	400mm	76 mm	≥100mm	71g
	450mm	84mm	≥115mm	71g
	500mm	92mm	≥135mm	71g